

NEW

HIGH POWER UV-VIS FIBER LIGHT SOURCE

L10290

OVERVIEW

The L10290 is a newly developed UV-visible fiber light source that outputs 200 nm to 1100 nm light through a light guide (sold separately).

Using a "high brightness lamp", the L10290 provides twice radiant intensity that of conventional models. Besides a compact size and light weight that make it simple to carry, the L10290 is specifically designed for easy use. These features make the L10290 ideal for all types of portable devices.




* Light guide is sold separately.

FEATURES

- **High Output: Twice Intensity**
(Compared to Conventional Model)
- **High Stability: Fluctuation 0.004 %p-p Typ.**
(Equivalent to 2×10^{-5} A.U.)
- **Long Life Lamp: 2000 Hours**
- **External Control**
- **Shutter Function**
- **Filter Holder (Sold Separately)**


SPECTRAL DISTRIBUTION (TYPICAL DATA)

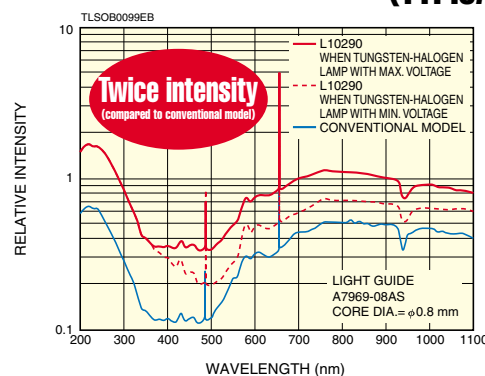

APPLICATIONS

- Spectrophotometry
- Environmental Measurement
- Absorption Spectrum Measurement
- Micro Titer Plate Reader
- Pharmaceutical Testing
- Film Thickness Measurement
- Semiconductor Inspection
- Optical Component Inspection

LIGHT OUTPUT STABILITY (at 230 nm)

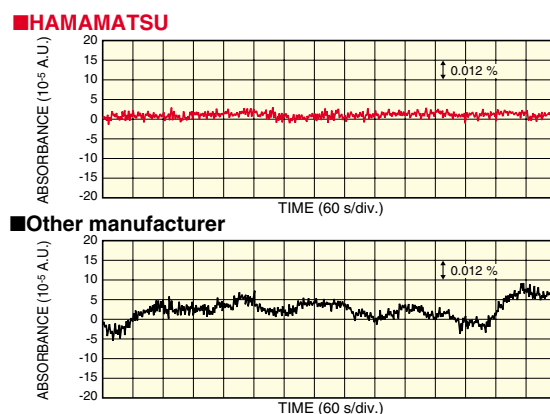

TLSOB0059EB

HAMAMATSU

SPECIFICATIONS

Parameter	Specifications
Spectral Distribution	200 nm to 1100 nm ①
Built-in Lamp	High Brightness Deuterium Lamp 9 W Tungsten Halogen Lamp ②
Output Stability	Fluctuation (Typ.) 0.004 % p-p (Equivalent to 2×10^{-5} A.U.) Drift (Max.) ± 0.3 %/h
Input Voltage (AC)	100 V to 240 V (100 V/200 V Auto Switching) Single Phase 50 Hz/60 Hz
Power Consumption (Max.)	140 VA
Weight	Approx. 4 kg
Light Guide Connector	SMA905D
Shutter	Included
External Control	Lamp ON/OFF ③
	Shutter open/close
	Status Signal for Lamp ON/OFF and Main Power ON/OFF
Main Unit Switch	Main Power ON/OFF
	Lamp ON/OFF ③
	Shutter open/close


NOTE: ①The deuterium lamp (200 nm to 400 nm) and tungsten halogen lamp (400 nm to 1100 nm) can each radiate light separately.
②Tungsten halogen lamp intensity is adjustable within a range of ± 20 %.
③The deuterium lamp and tungsten halogen lamp can each be turned on or off separately.

REPLACEMENT LAMP

Type No.	Lamp	Lamp Life ④
L10295	High Brightness Deuterium Lamp	2000 h
L10296	9 W Tungsten Halogen Lamp	2000 h

NOTE: ④Guaranteed life. Defined as the time when the radiant intensity falls to 50 % of its initial value.

DIMENSIONAL OUTLINES (Unit: mm)

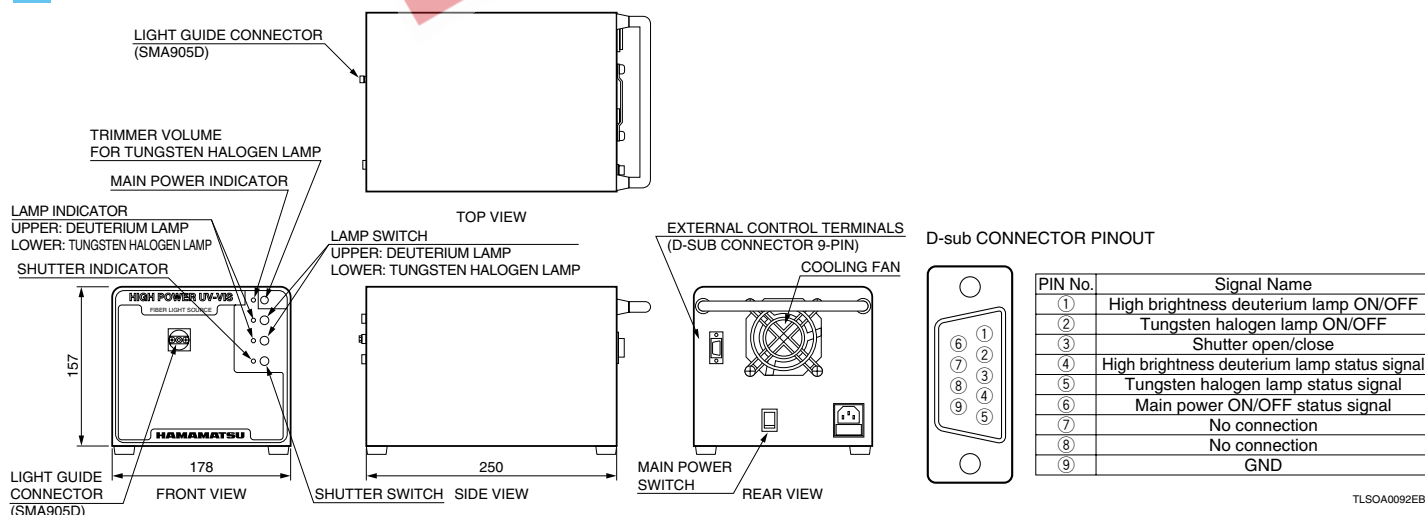

Subject to local technical requirements and regulations, availability of products included in this promotional material may vary. Please consult with our sales office. Information furnished by HAMAMATSU is believed to be reliable. However, no responsibility is assumed for possible inaccuracies or omissions. Specifications are subject to change without notice. No patent rights are granted to any of the circuits described herein. ©2006 Hamamatsu Photonics K.K.


LIGHT GUIDE ⑤ (SOLD SEPARATELY)

Type No.	Core Diameter	Length
A7969-06AS ⑥	0.6 mm	1.0 m
A7969-08AS ⑥	0.8 mm	

NOTE: ⑤Other shapes and combinations are available upon request.

⑥AS type is a UV resistant light guide.

UV RESISTANT CHARACTERISTICS (at 215 nm)

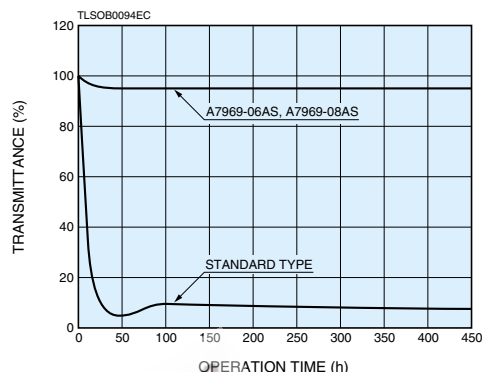

FILTER HOLDER (SOLD SEPARATELY)

Type No.	Applicable Filter
A10394	3 mm thick, 1 inch dia. or 1 × 1 inch

Filters should be prepared by the user. When attaching a filter holder, remove the cover and insert the filter holder into the lamp house.

HAMAMATSU

WEB SITE www.hamamatsu.com

HAMAMATSU PHOTONICS K.K., Electron Tube Division

314-5, Shimokanzo, Iwata City, Shizuoka Pref., 438-0193, Japan, Telephone: (81)539/62-5248, Fax: (81)539/62-2205

U.S.A.: Hamamatsu Corporation, 360 Foothill Road, P. O. Box 6910, Bridgewater, N.J. 08807-0910, U.S.A., Telephone: (1)908-231-0960, Fax: (1)908-231-1218 E-mail: usa@hamamatsu.com

Germany: Hamamatsu Photonics Deutschland GmbH: Arzbergerstr. 10, D-82211 Herrsching am Ammersee, Germany, Telephone: (49)8152-375-0, Fax: (49)8152-2658 E-mail: info@hamamatsu.de

France: Hamamatsu Photonics France S.A.R.L.: 19, Rue du Saule Trappu, Parc du Moulin de Massy, 91882 Massy Cedex, France, Telephone: (33)1 69 53 71 00, Fax: (33)1 69 53 71 10 E-mail: infos@hamamatsu.fr

United Kingdom: Hamamatsu Photonics UK Limited: 2 Howard Court, 10 Tewin Road Welwyn Garden City Hertfordshire AL7 1BW, United Kingdom, Telephone: 44-(0)1707-294888, Fax: 44(0)1707-325777 E-mail: info@hamamatsu.co.uk

North Europe: Hamamatsu Photonics Norden AB: Smidesvägen 12, SE-171-41 SOLNA, Sweden, Telephone: (46)8-509-031-00, Fax: (46)8-509-031-01 E-mail: info@hamamatsu.se

Italy: Hamamatsu Photonics Italia: S.R.L.: Strada della Moia, 1/E, 20020 Arese, (Milano), Italy, Telephone: (39)02-935 81 733, Fax: (39)02-935 81 741 E-mail: info@hamamatsu.it

TL501059E02
OCT. 2006 IP